

Remove unwanted pigment
and enjoy smoother skin!

- Benign Pigmentation Treatment
- Various Color Tattoo Removal
- Scar Treatment
- Skin Rejuvenation

PICOCARE 450



Before & After

Melasma



| Before |



| After |

After 2 Tx Time
L532 nm, 4 mm, 0.3 J/cm², 2 Hz

ABNOM



| Before |



| After |

After 2 Tx Time
1064 nm, 10 mm, 0.8 J/cm², 10 Hz

Body Tattoo




| Before |



| After |

After 2 Tx Time
1064 nm, 3 mm, 2.5 J/cm², 2 Hz

Large Pores

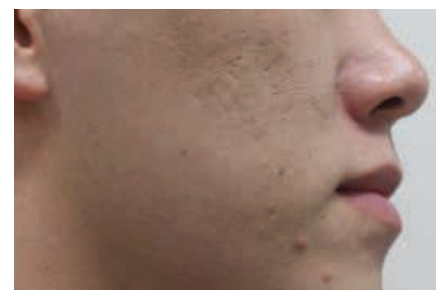

| Before |



| After |

After 1 Tx Time
1064 nm HEXA, 10 mm, 0.3 J/cm²,
10 Hz

PICOCARE 450

Picosecond Laser

Irradiation type	Nd:YAG
Wavelength	1064 nm & 532 nm + 595, 660 nm
Max. energy	600 mJ (1064 nm), 200 mJ (532 nm)
Pulse duration	450ps
Peak power	~ 1.33 GW (1064 nm), ~ 0.44 GW (532 nm)
Spot size	2 ~ 10 mm
Repetition rate	Single, 1 ~ 10 Hz
Delivery	Articulated Arm
Optional handpiece	Dye H/P : 595, 660 nm
Dimension	450 x 939 x 908 (W x D x H) mm
Weight	90 kg / 4kVA



Head Office
64 Techno 8-ro, Yuseong-gu, Daejeon, Korea

WON TECH INC
500 W Office Center Drive Suite #400-4108 Fort Washington,
PA 19034, USA
(+1) 267-342-8890 james@wtlaser.com

WON TECH JAPAN CO., LTD.
1F, 3-16-6, Hatchobori, Chuo-ku, Tokyo, Japan 104-0032
(+81) 3-5542-8566 nhan@wtlaser.com

WON TECH CHINA CO., LTD.
F510, Block C, Tower 1, Wangjing Soho, Wangjing Avenue,
Chaoyang Dist., Beijing, 100102, P.R. CHINA
(+86) 010-6464-3241 bjwtk@wtlaser.com.cn

Pangyo Branch Office
3F, 22, Deawangpangyo-ro 712beon-gil, Bundang-gu,
Seongnam-si, Gyeonggi-do, Korea

E-mail. globalsales@wtlaser.com
Tel (+82) 31-1670-1450 Fax (+82) 42-934-9491

www.wtlaser.com

f www.facebook.com/wtlaser

www.instagram.com/wontech_laser



For Export Only

PICOCARE 450

Picosecond Nd:YAG laser
1064nm, 532nm,
595nm & 660nm

Beyond Nano,
Experience Real PICO





What is PICOCARE 450?

Innovative Technology for Optimal Pigmentation Treatment

PICOCARE 450 delivers better photomechanical effect on the targeting pigmentation or unwanted tattoo particles with less thermal effects on the surrounding skin tissues. This is because of higher peak power of ultra-short pulses and picoseconds that are 1000 times shorter (one trillionth of a second) than conventional Q-switched lasers.


Scar Treatment Available by PICOCARE 450

PICOCARE 450 with HEXA MLA can effectively treat atrophic scars and hypertrophic scars with less discomfort and low rate of side effects, replacing aggressive scar treatments such as CO₂ laser and subcision.


NANO VS PICO



[Conventional Nanosecond Laser]


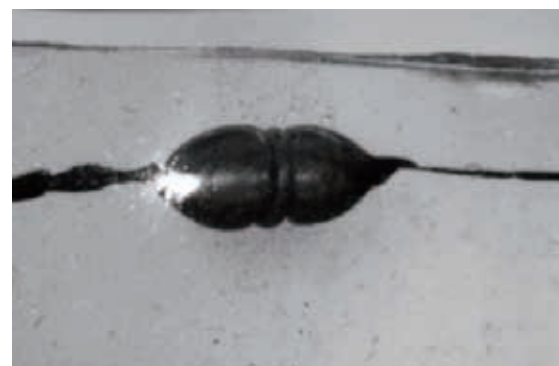


[PICOCARE 450]

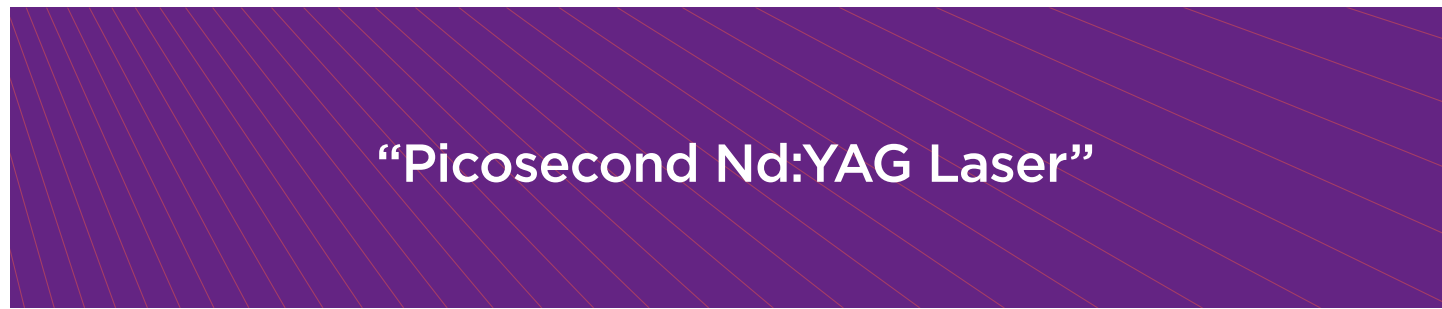


Multiple Wavelengths

PICOCARE 450 uses 1064 nm, 532 nm, 595 nm & 660 nm wavelengths to remove multiple benign pigmentations of the epidermis and dermis, and unwanted various color tattoo.

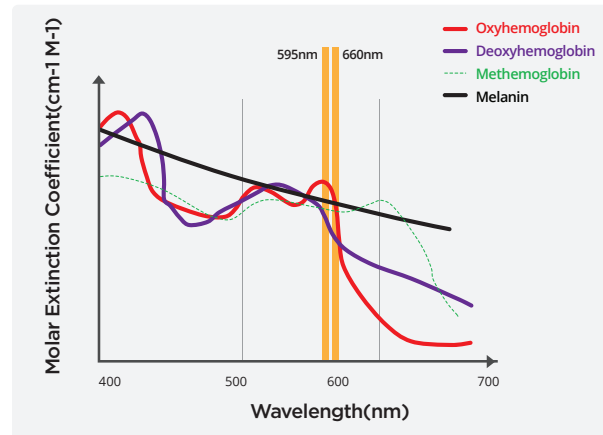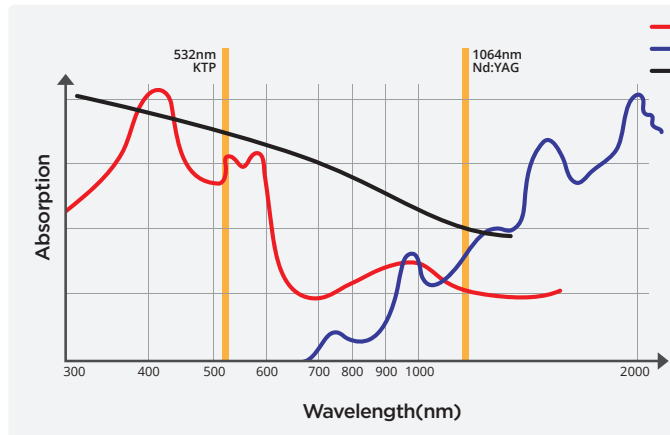


Picocare 450 Handpiece

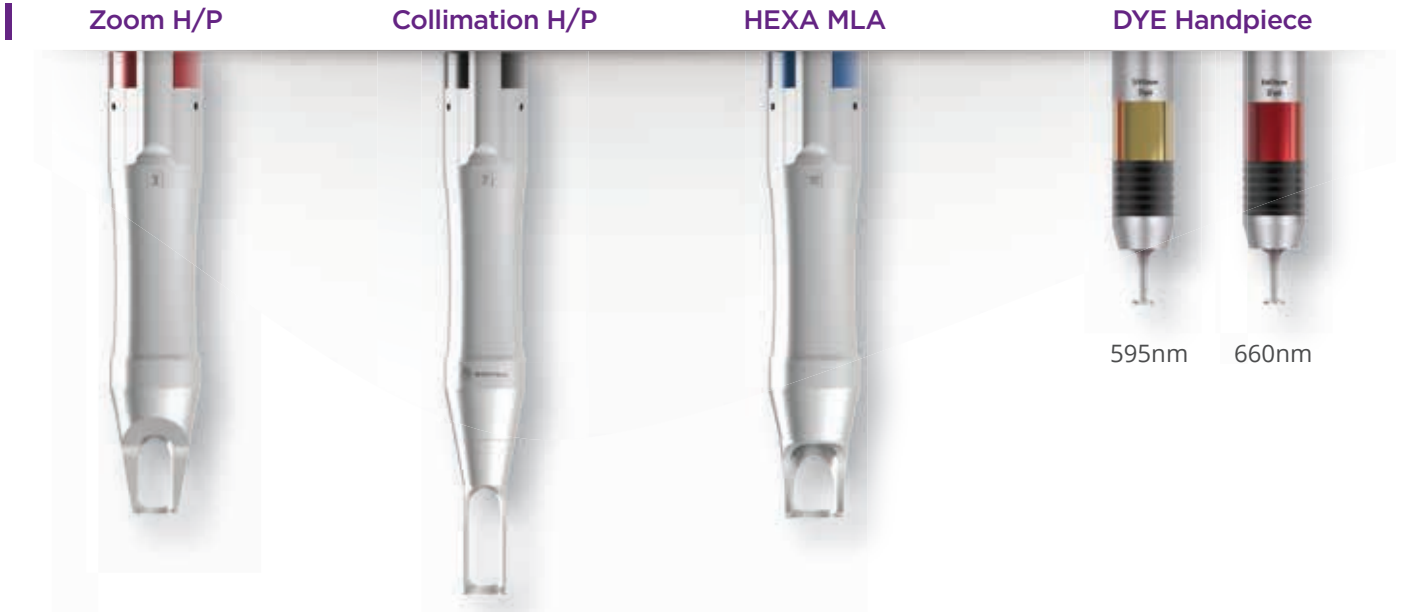

HEXA MLA

PICOCARE 450 can also treat scars and skin rejuvenation by using HEXA MLA. The mechanism is that LIOB (Laser Induced Optical Breakdown) is created on the target area, stimulating collagen and elastic fibers with stratum corneum preserved.

Only 1 session of scar treatment can show dramatic results in volumized atrophic scars with fast healing process and relatively short down time.

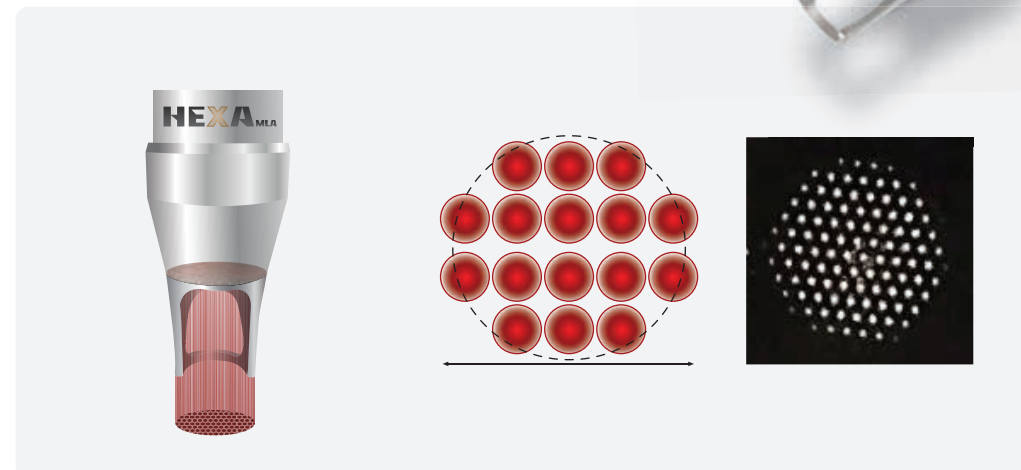

Smart and Convenient Laser System

- The compact size compared to other picosecond lasers
- Real-time auto calibration
- Fast initiation time within 2 minutes
- Intuitive user interface

Indications

Atrophic Scars & Hypertrophic Scars	Skin Rejuvenation	Epidermal Benign Pigmented Lesions	Dermal Benign Pigmented Lesions	Various Color Tattoo
<ul style="list-style-type: none"> - Acne Scars - Chicken Pox - Nail Scratches - Stretch Marks 	<ul style="list-style-type: none"> - Uneven Skin Tone - Dull Skin Texture - Large Pores - Fine Wrinkles 	<ul style="list-style-type: none"> - Freckles - Café-au-lait - Lentigines - Age spot 	<ul style="list-style-type: none"> - ABNOM - Ota nevus - Ito Ota Blue - Melasma 	<ul style="list-style-type: none"> - Black - Blue - Green - Yellow - Red